


pDendra2-C vector

The vector sequence has been compiled using the information from sequence databases, published literature, and other sources, together with partial sequences obtained by Evrogen. This vector has not been completely sequenced.

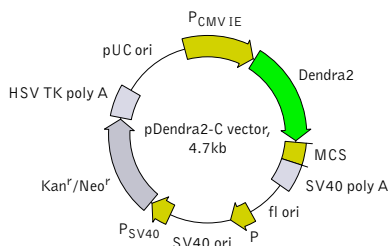

For vector sequence, please visit our Web site at <http://www.evrogen.com/support/vector-info.shtml>

pDendra2-C vector MCS

Dendra2
 ... TCC. GGC. GAC. AGC. GGC. GTG. TAC. AAG. ACT. CGA. GCT. CAA. GCT. TCG. AAT. TCT. GCA. GTC. GAC. GGT. ACC. GCG. GGC. CCG. GGA. TCC. ACC. GGA. TCT. AGA. TAA. CTG. ATC. A. . .
 * — not unique sites.
 # — sites are blocked by *dam* methylation. If you wish to digest the vector with these enzymes, you will need to transform the vector into a *dam*⁻ host and make fresh DNA.

Location of features

P_{CMV IE}: 1-589
 Enhancer region: 59-465
 TATA box: 554-560
 Transcription start point: 583
 Dendra2
 Kozak consensus translation initiation site: 606-616
 Start codon (ATG): 613-615; Stop codon: 1393-1395
 Last amino acid in Dendra: 1300-1302
 MCS: 1327-1393
 SV40 early mRNA polyadenylation signal
 Polyadenylation signals: 1535-1540 & 1564-1569
 mRNA 3' ends: 1573 & 1585
 f1 single-strand DNA origin: 1632-2087
 Eukaryotic promoter for expression of Kan^r gene
 -35 region: 2149-2154; -10 region: 2172-2177
 Transcription start point: 2184
 SV40 origin of replication: 2428-2563
 SV40 early promoter
 Enhancer (72-bp tandem repeats): 2261-2332 & 2333-2404
 21-bp repeats: 2408-2428, 2429-2449 & 2451-2471
 Early promoter element: 2484-2490
 Major transcription start points: 2480, 2518, 2524 & 2529
 Kanamycin/neomycin resistance gene
 Neomycin phosphotransferase coding sequences:
 Start codon (ATG): 2612-2614; Stop codon: 3404-3406
 G->A mutation to remove Pst I site: 2794
 C->A (Arg to Ser) mutation to remove BssH II site: 3140
 Herpes simplex virus (HSV) thymidine kinase (TK) polyadenylation signal
 Polyadenylation signals: 3642-3647 & 3655-3660
 pUC plasmid replication origin: 3991-4634

References

Gorman (1985). "High efficiency gene transfer into mammalian cells." In: *DNA cloning: A Practical Approach, Vol. II*. Ed. by Glover. (IRL Press, Oxford, U.K.) Pp. 143-90.

Haas et al. (1996) "Codon usage limitation in the expression of HIV-1 envelope glycoprotein." *Curr Biol*, 6(3): 315-24 / pmid: 8805248

Kozak (1987) "An analysis of 5'-noncoding sequences from 699 vertebrate messenger RNAs." *Nucleic Acids Res*, 15(20): 8125-48 / pmid: 3313277

Product	Cat.#	Size
pDendra2-C vector	FP821	20 µg
The price does not include delivery. The price varies in different countries. Please contact your local distributor for exact prices and delivery information.		
Vector type	mammalian expression vector	
Reporter	Dendra2	
Reporter codon usage	mammalian	
Promoter for Dendra2	P _{CMV IE}	
Host cells	mammalian	
Selection	prokaryotic - kanamycin eukaryotic - neomycin (G418)	
Replication	prokaryotic - pUC ori eukaryotic - SV40 ori	
Use	Dendra2 expression in mammalian cells; generation of fusions to the Dendra2 C-terminus	

Vector description

pDendra2-C is a mammalian expression vector encoding green/red fluorescent protein Dendra2. The vector allows generation of fusions to the Dendra2 C-terminus and expression of Dendra2 fusions or Dendra2 alone in eukaryotic (mammalian) cells.

Dendra2 codon usage is optimized for high expression in mammalian cells (humanized) [Haas et al. 1996]. To increase mRNA translation efficiency, Kozak consensus translation initiation site is generated upstream of the Dendra2 sequence [Kozak 1987]. Multiple cloning site (MCS) is located between Dendra2 coding sequence and SV40 polyadenylation signal (SV40 polyA).

The vector backbone contains immediate early promoter of cytomegalovirus (P_{CMV IE}) for protein expression, SV40 origin for replication in mammalian cells expressing SV40 T-antigen, pUC origin of replication for propagation in *E. coli* and f1 origin for single-stranded DNA production. SV40 polyadenylation signals (SV40 poly A) direct proper processing of the 3'-end of the reporter mRNA.

SV40 early promoter (P_{SV40}) provides neomycin resistance gene (Neo^r) expression to select stably transfected eukaryotic cells using G418. Bacterial promoter (P) provides kanamycin resistance gene expression (Kan^r) in *E. coli*. Kan^r/Neo^r gene is linked with herpes simplex virus (HSV) thymidine kinase (TK) polyadenylation signals.

Generation of Dendra2-fusion proteins

A localization signal (or a gene of interest) should be cloned into MCS of the vector. It will be expressed as a fusion to the Dendra2 C-terminus when inserted in the same reading frame as Dendra2 and no intervening stop codons are present. Dendra2-tagged fusions retain fluorescent properties of the native protein allowing fusion localization *in vivo*. Unmodified vector will express Dendra2, when transfected into eukaryotic (mammalian) cells.

Note: The plasmid DNA was isolated from *dam*⁺-methylated *E. coli*. Therefore some restriction sites are blocked by methylation. If you wish to digest the vector using such sites you will need to transform the vector into a *dam*⁻ host and make fresh DNA.

Expression in mammalian cells

pDendra2-C vector can be transfected into mammalian cells by any known transfection method. If required, stable transformants can be selected using G418 [Gorman 1985].

Propagation in *E. coli*

Suitable host strains for propagation in *E. coli* include DH5alpha, HB101, XL1-Blue, and other general purpose strains. Plasmid incompatibility group is pMB1/ColE1. The vector confers resistance to kanamycin (30 µg/ml) to *E. coli* hosts. Copy number in *E. coli* is about 500.

Notice to Purchaser:

Photoactivatable FP-related products are intended to be used by academic (non-commercial) entities and for research purposes only. Any use of the proprietary nucleic acid or protein other than for research use or by a commercial entity is strictly prohibited. Transfer of this product by purchaser to any other party is specifically prohibited.

The CMV promoter is covered under U.S. Patents 5,168,062 and 5,385,839, and its use is permitted for research purposes only. Any other use of the CMV promoter requires a license from the University of Iowa Research Foundation, 214 Technology Innovation Center, Iowa City, IA 52242.

MATERIAL SAFETY DATA SHEET INFORMATION: To the best of our knowledge, these products do not require a Material Safety Data Sheet. However, all the properties of these products (and, if applicable, each of their components) have not been thoroughly investigated. Therefore, we recommend that you use gloves and eye protection, and wear a laboratory coat when working with these products.